

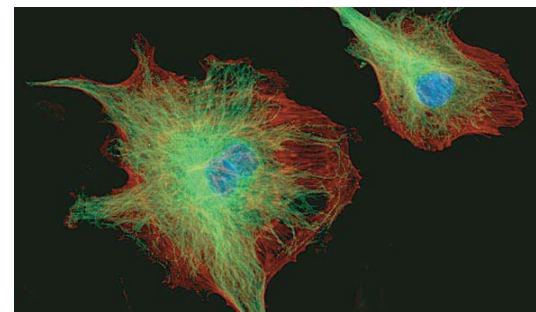


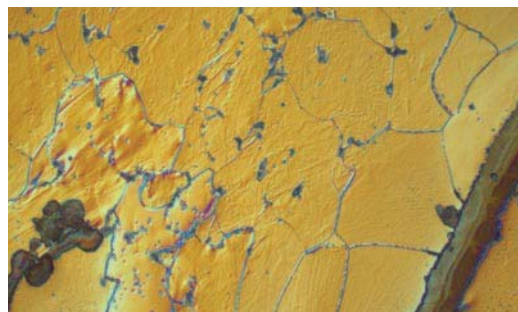

## ProgRes® CapturePro Software

Advance your image results

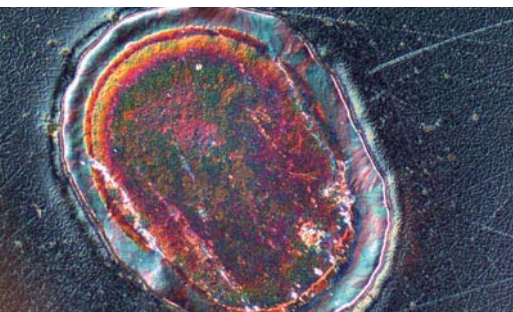

### Strong software - simply ingenious

Whether for simple image documentation or demanding image analysis: with merely three mouse clicks you can optimize your images - thanks to the automatic functions included in ProgRes® CapturePro software.

#### Benefits

- CapturePro is included in the delivery of ProgRes® microscope cameras
- For Microsoft Windows and Apple Macintosh OS
- Free updates for registered users
- Stand-alone software and TWAIN Plug-In
- Optional software development kit (SDK)

[ FULL WIN 7 Support ]  
32 and 64 bit

### Optimal image quality for ProgRes® cameras

CapturePro supports the excellent quality of ProgRes® microscope cameras and delivers optimal image quality and reproducible results. Calibrations and image correction matrices, for example, shadings can be individually generated for your camera and are saved in your own user profile.


#### Features

- Full support of *SpeedXT<sup>core</sup>* & Multicore-processors
- Automatic white balance
- Multi-fluorescence mode
- Measurement capabilities & Image annotation
- Multi-focus
- Time-lapse
- Continuous automatic exposure control
- Shading calibrations
- Multi-camera operation
- Individual user profiles


# ProgRes® CapturePro Software

## Advance your image results

### Features & Specifications

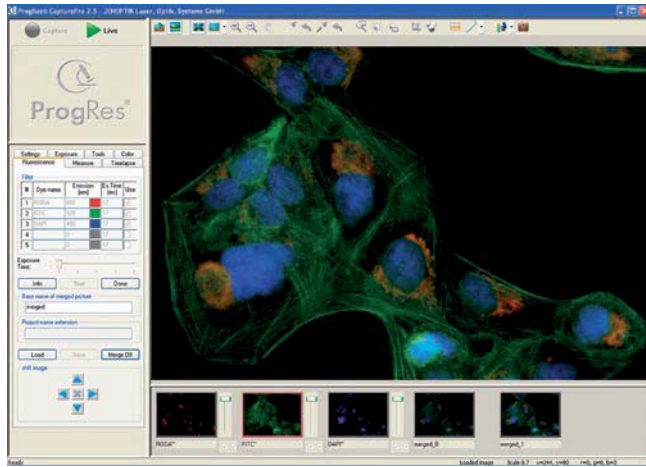

Fluorescence mode: merging of multi-fluorescence images

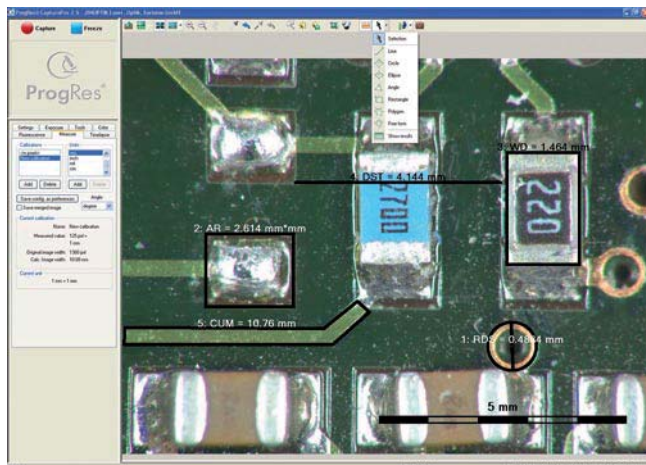

Easy-to-use tools for measurement and annotation

#### Multicore-Processor support

- Jenoptik image processing pipe is designed to use the fully available core processor power

#### Shading calibration

- Noise can be reduced by black shading calibration
- Vignetting of lenses can be eliminated by white shading calibration

#### User Preferences and Profiles

- Storage of individual user settings
- Simultaneous connection of up to 5 cameras to one computer

#### Multi-focus (Z-Stacking)

- Capture images from several planes within live image
- Multi-focus combines the sharp areas of each plane to create an image that is sharp throughout

#### Time lapse

- Adjustment of delay and duration of capture sequences
- Storage as movie file

#### Continuous automatic exposure

- The ideal exposure time will be calculated integrally over the entire image automatically

#### Multi-Fluorescence mode

- Definition of up to 5 color filters
- Guided workflow for capturing multi-fluorescence images
- DIC-images overlay

#### Measurements

- Performed in live-image mode
- Calibration options
- Various measurement tools
- Annotations

#### Live-image functionality


- Rotate/ mirror
- Focus tool
- Automatic exposure control
- Automatic white balance can be activated with one button click over the full live image

#### Further features

- Text stamp
- Full screen mode
- Print button
- Snap shot
- Programmable resolution

#### Hardware requirements

- PC: MS WIN XP / Win Vista / WIN 7
- 32 and 64 Bit Systems
- 3 GHz CPU, 1 GB RAM, 64 MB graphics, FireWire a or USB 2.0, Multicore recommended

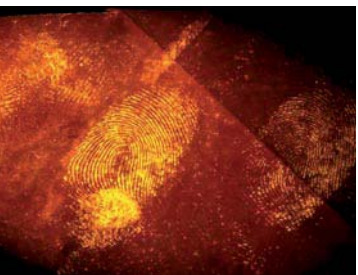

It is our policy to constantly improve the design and specifications of our products. Accordingly, the details represented herein cannot be regarded as final and binding.



**Ryf AG**  
 Bettlachstrasse 2  
 2540 Grenchen  
 tel 032 654 21 00  
 fax 032 654 21 09  
[www.ryfag.ch](http://www.ryfag.ch)