


Mx Link

Data export through WIFI

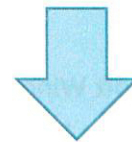


MxLink is a Microsoft Windows application for the data export through WIFI from the Mx200 to a PC or Network drive.

The MxLink application allows users to configure one or more Mx200 devices to send measurement results directly to a specified folder on the target PC. This folder can be local to the MxLink PC, or a mapped network folder, accessible by the MxLink PC.

Data received in the MxLink application will be displayed in a separate "tabbed" window for each device that is connected. This data of the last export can be sent directly to a Microsoft Excel sheet, or it can be copied to, and pasted from, the Windows clipboard.

| Name | Coef | Nominal | Actual | Tol- | Tol+ | Deviation | Tendency |
|------------|------|---------|--------|---------|--------|-----------|----------|
| Point 1 | X | 2,100 | 2,140 | | | | |
| Point 1 | Y | 1,600 | 1,600 | | | | |
| Circle 2 | X | 5,223 | 5,223 | | | | |
| Circle 2 | Y | 7,148 | 7,148 | | | | |
| Circle 2 | D | 5,600 | 5,633 | -0,100 | 0,100 | 0,033 | 33 |
| Line 3 | X | 4,200 | 4,200 | | | | |
| Line 3 | Y | 1,110 | 1,110 | | | | |
| Line 3 | B | 0,100° | 0,000° | -1,000° | 1,000° | -0,100° | -10 |
| Distance 4 | X | 0,000 | 0,000 | | | | |
| Distance 4 | Y | 6,000 | 6,038 | -0,100 | 0,100 | 0,038 | 38 |
| Distance 4 | L | 6,038 | 6,038 | | | | |

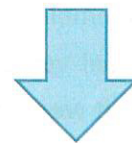

| Name | Coef | Nominal | Actual | Tol- | Tol+ | Deviation | Tendency |
|------------|------|---------|--------|---------|--------|-----------|----------|
| Point 1 | X | 2,1 | 2,14 | | | | |
| Point 1 | Y | 1,6 | 1,6 | | | | |
| Circle 2 | X | 5,223 | 5,223 | | | | |
| Circle 2 | Y | 7,148 | 7,148 | | | | |
| Circle 2 | D | 5,6 | 5,633 | -0,1 | 0,1 | 0,033 | 33 |
| Line 3 | X | 4,2 | 4,2 | | | | |
| Line 3 | Y | 1,11 | 1,11 | | | | |
| Line 3 | B | 0,100° | 0,000° | -1,000° | 1,000° | -0,100° | -10 |
| Distance 4 | X | 0 | 0 | | | | |
| Distance 4 | Y | 6 | 6,038 | -0,1 | 0,1 | 0,038 | 38 |
| Distance 4 | L | 6,038 | 6,038 | | | | |

MICROSCOPY + METROLOGY SERVICES
Suisse made.



RYF AG
Bettlachstrasse 2
CH-2540 Grenchen
Tel. +41 32 654 21 00
ryfag@ryfag.ch · www.ryfag.ch

RYF SA Succursale
Route de Genève 9c
CH-1291 Commugny
Tel +41 22 776 82 28
ryfag@ryfag.ch · www.ryfag.ch

RYF AG (Zürich)
Bahnhofplatz 17
CH-8400 Winterthur
Tel +41 52 560 22 25
ryfag@ryfag.ch · www.ryfag.ch